

## Hydrogen-Deuterium Spectrum Experiment

### 1. Introduction

In this experiment, we used the wavelength differences between hydrogen and deuterium Balmer series to calculate the mass ratio of hydrogen and deuterium.

First, we use a monochromator to scan over the sodium spectrum at a constant rate. Since the wavelength difference between the two lines of the yellow sodium doublet is known (sodium D lines), it can be used to calculate a conversion factor between wavelength difference and time difference that depends on the scanning rate we choose. Then we measured the wavelength differences between hydrogen and deuterium Balmer series (namely:  $\alpha$ ,  $\beta$ ,  $\gamma$ ,  $\delta$ , and  $\epsilon$  lines). Using the expression of effective mass and orbital energy difference

$$\mu = \frac{mM}{M+m}, \quad E_i - E_f \propto \mu \left( \frac{1}{n_f^2} - \frac{1}{n_i^2} \right)$$


we find that

$$\frac{\Delta\lambda_{HD}}{\lambda_H} = \frac{\lambda_H - \lambda_D}{\lambda_H} = \frac{1/\mu_H - 1/\mu_D}{1/\mu_H} = \frac{1 - M_H/M_D}{1 + M_H/m}$$

Since  $M_H/m$  and  $\lambda_H$  is known and we can measure  $\Delta\lambda_{HD}$ , we then can solve for the desired quantity  $M_H/M_D$ . To improve accuracy, we plot  $\Delta\lambda_{HD}$  for all lines in the Balmer series, and fit the data to find the slope (i.e.  $\frac{\Delta\lambda_{HD}}{\lambda_H}$ ) instead of just using one single spectral line.

### 2. Result for the mass ratio

The data plot for  $\Delta\lambda_{HD}$  with respect to hydrogen wavelength is shown below (note that the fitted line is included also)

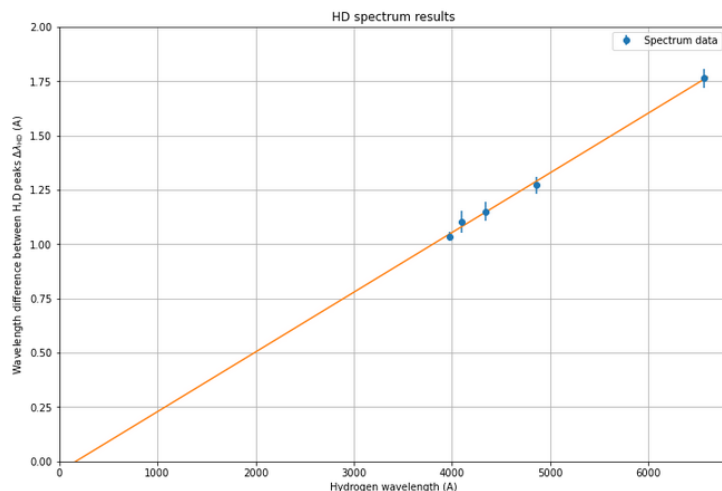

The fitted line gives a slope of  $0.000257 \pm 0.000009$ . Using the equation given above, the mass ratio is:

$$\frac{M_H}{M_D} = 0.50 \pm 0.02$$

The accepted value of hydrogen-deuterium mass ratio is  $M_H/M_D = 0.500248$ , which is only about 0.0002 off from our measured value, so lies within the uncertainty (0.02). This concluded that our measured value agree very well with the accepted value.

### 3. Possible error

Assume we made a mistake in measurements by pausing the wavelength scanning for 2 seconds between two sodium D lines, how would it effects our mass ratio final result?

First, we consider how this mistake will effect the scaling factor  $K_{Na} = \frac{\Delta\lambda_D}{\Delta t_D}$ . Since the time interval between sodium peaks is lengthened by our mistake

$$K_{Na} = \frac{\Delta\lambda_D}{\Delta t_D} \rightarrow \frac{\Delta\lambda_D}{\Delta t_D + 2s}$$

The scaling factor will become smaller, namely, it is scaled by  $\frac{\Delta t_D}{\Delta t_D + 2s} = 0.94026 \pm 0.00004$  (with  $\Delta t_D = 31.476s \pm 0.025$  from our analysis).

Since we calculate  $\Delta\lambda_{HD}$  values using  $K_{Na}$  as follow

$$\Delta\lambda_{HD} = K_{Na}\Delta t_{HD}$$

the new  $\Delta\lambda_{HD}$  will also be scaled by the same amount (i.e  $\Delta\lambda_{HD} \rightarrow 0.94026\Delta\lambda_{HD}$ ). So the slope ( $\frac{\Delta\lambda_{HD}}{\lambda_H}$ ) will also decrease by a factor of 0.94026. Solving for the mass ratio in the equation on the first page, with the known value  $M_H/m = 1836.15$ , we get

$$\frac{M_H}{M_D} = 1 - 1837.15 \times \frac{\Delta\lambda_{HD}}{\lambda_H}$$

this means that a 0.94026 decrease in slope will result in a new mass ratio of  $0.525 \pm 0.015$  (original slope is  $0.000275 \pm 0.000009$ , detailed calculation is done in the analysis notebook). So, by taking the ratio of the "correct mass ratio" and the "incorrect mass ratio"

$$\frac{(M_H/M_D)_{incorrect}}{(M_H/M_D)_{correct}} = \frac{0.525 \pm 0.015}{0.50 \pm 0.02} = 1.05 \pm 0.0516$$

So the mistake we made causes us to measure a mass ratio that is 105% of the mass ratio without making the mistake.